PharmaJet[®] Tropis[®]

EN-Needle-free Injection System

DE-Nadelloses Injektionssystem

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Volume : 0,1 ml (+/-10 %)

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Symbols Glossary

(Note: All referenced symbols are derived from **ISO 15223-1**, *Medical Devices - Symbols to be used with medical device labels, labelling and information to be supplied.*)

Symbol	Title of Symbol	Description of Symbol	Ref. No.
	Manufacturer	Indicates the medical device manufacturer.	5.1.1
	Authorized representative in the European Community	Indicates the authorized representative in the European Community.	5.1.2
	Use-by date	Indicates the date after which the medical device is not to be used.	5.1.4
	Batch code	Indicates the manufacturer's batch code so that the batch or lot can be identified.	5.1.5
	Catalogue number	Indicates the manufacturer's catalogue number so the batch or lot can be identified.	5.1.6
	Sterilized using ethylene oxide	Indicates a medical device that has been sterilized using ethylene oxide.	5.2.3
	Do not use if package is damaged	Indicates a medical device that should not be used if the package has been damaged or opened.	5.2.8
	Temperature limit	Indicates the temperature limits to which the medical device can be safely exposed.	5.3.7
	Do not re-use	Indicates a medical device that is intended for one use, or for use on a single patient during a single procedure.	5.4.2
	Consult instructions for use	Indicates the need for the user to consult the instructions for use.	5.4.3
	Caution	Indicates the need for the user to consult the instructions for use for important cautionary information such as warnings and precautions that cannot, for a variety of reasons, be presented on the medical device itself.	5.4.4
	Non-pyrogenic	Indicates a medical device that is non-pyrogenic.	5.6.3
	Do not touch	Indicates components that should not be touched.	N/A

PharmaJet[®]

Tropis[®]

English

System Components

Injector

Claimed Lifetime: 20,000 cycles

Syringe

STERILE EO

Filling Adapter

STERILE EO

CAUTION: Federal law restricts this device to sale by or on the order of a physician.

Indications for Use

The PharmaJet Needle-free Injection System is intended to deliver various medications and vaccines intradermally by means of a narrow, high velocity fluid jet, which penetrates the skin and delivers the medicine or vaccine into the dermis. Healthcare providers who routinely administer injections may use the PharmaJet Needle-free Injection System. It may also be used by patients authorized by their physicians to self-inject, or have other individuals administer injections of prescribed medication.

Clinical Notes

The following local responses may occur with Needle or Needle-free Injection Technology.

- Some medications may produce a burning, aching, or stinging sensation. Check the medication information for a description of local reactions.
- Following the injection, a drop of blood may appear. Blood and/or fluid loss can be minimized by immediately applying pressure to the injection site for approximately 60 seconds with a cotton ball or gauze. Apply adhesive bandage if needed or desired.
- To ensure delivered dose to target tissue, refer to **Section 4) Give Injection** for proper technique. For a partial or incomplete injection, refer to the medication or vaccine product labeling.
- In the event of changes in performance of the device, DO NOT USE, contact PharmaJet.

Cleaning, Disinfecting, Storage, and Disposal

CAUTION: DO NOT autoclave.


CLEANING & DISINFECTING: The Injector may be cleaned by wiping the exposed surfaces using 70% Isopropyl Alcohol.

STORAGE: Ensure that the Hammer is in the uncharged position.

NOTE: The Injector is resistant to short term duration exposure to water.

DISPOSAL: Dispose of used devices in accordance with local requirements.

**Storage and Transport
Temperature:**

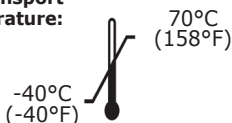
Reorder Information

<u>Item</u>	<u>Catalog #</u>
Injector	C503-10
Syringes – 100 Single Units	C503-34
Syringes – 100 Bulk Units	C503-35
Syringes – 200 Bulk Units	C503-36
Filling Adapters – 100 Single Units	C503-42

CAUTION: Use only PharmaJet components with the Tropis® Needle-free Injection System.

**Operating
Temperature:**

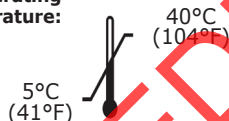Warranty Information

The Injector carries a replacement warranty against defects in materials and workmanship for two years from the date of purchase.

Contact Information

PharmaJet®

US Office:

PharmaJet, Inc.
400 Corporate Circle, Suite N
Golden, Colorado 80401
USA

Tel: +1 888 900 4321
Tel Outside the U.S.: +1 303 526 4278

EC REP

EU Authorized Representative:

Emergo Europe
Prinsessegracht 20
2514 AP The Hague
The Netherlands

Warnings and Cautions

1) Prepare Injector

 **CAUTION:** Repeated dry firing will damage Injector. Dry firing is activating the Injector with either an empty Syringe or without a Syringe.

 **CAUTION:** DO NOT open Injector Handles when device is charged.

2) Fill Syringe

 **WARNING:** The tip of the Filling Adapter is sharp. Use caution when handling.

 **WARNING:** Syringes auto-disable during use. Attempting to use a Syringe more than one time may result in contamination or device malfunction.

 **WARNING:** Vial integrity or sterility is not ensured when used with the Filling Adapter. DO NOT reuse the Filling Adapter in another vial. Doing so may contaminate contents of new vial.

 **CAUTION:** The protective Syringe Cap is designed to protect against touch contamination only, and does not provide a barrier against other forms of contamination.

 **CAUTION:** DO NOT attempt to use a partially filled Syringe; damage may occur.

 **CAUTION:** Visually inspect sterile package for damage. DO NOT use if torn or damaged. DO NOT re-sterilize.

3) Load Injector

 **WARNING:** DO NOT direct the Syringe toward the face or nearby skin when loading the Syringe into an Injector, preparing to give an injection, or ejecting a Syringe after use. Injury may occur.

 **CAUTION:** Partial insertion of a Syringe into an Injector may damage the Syringe and/or Injector.

 **CAUTION:** DO NOT load a Syringe if the Inner Core Housing is not visible and locked in position. This may cause inadvertent activation.

4) Give Injection

 **WARNING:** Improper injection site location or technique may result in injury.

 **CAUTION:** If the Injector fails to reset or if the retractable core does not move freely, DO NOT USE. Continued use may result in an incomplete injection or failure to inject.

 **CAUTION:** Patients receiving anticoagulants or persons with hemophilia or other coagulation disorders may have more post-injection bleeding and/or bruising.

 **CAUTION:** DO NOT deliver injection with the Injector at an angle other than 90 degrees to the patient, as shown in figure 11.

Suitable Injection Sites

 **WARNING:** DO NOT position the Syringe directly over a blood vessel. Although injection into the vessel is highly unlikely, injection over a vessel may produce an adverse reaction including, but not limited to, bleeding, bruising, swelling, or tenderness.


Select a site that is free from irritation, with a low likelihood of infection, and where the skin is intact and free from abrasions, bruises, etc.

Needle-free injection sites are identical to those used for conventional needle and syringe; the deltoid, suprascapular, and lateral thigh areas are most often used as the delivery sites for vaccine. Other types of injectant may have specific sites for injection. Follow package directions from the manufacturer of injectant.

Component Description


Injector

1. Activation Button
2. Inner Core Housing
3. Hammer
4. Eject Button
5. Safety Button
6. Handles
7. Latch Button

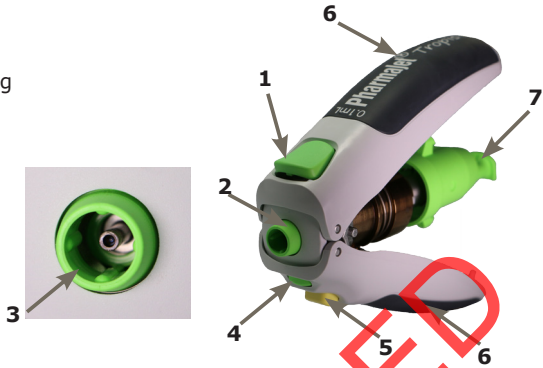
Syringe

1. Syringe Cap
2. Syringe Barrel
3. Syringe Plunger
4. Fill Line

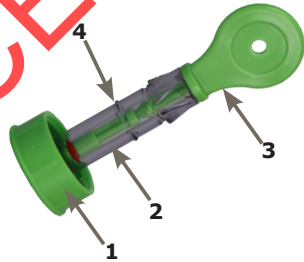Filling Adapter

1. Syringe Interface
2. Filling Adapter Spike

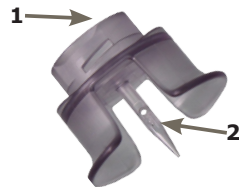Workflow

1) Prepare Injector

- a) Verify if the Injector is in the charged or uncharged state by confirming the position of the Injector Hammer as shown in figures 1a (uncharged) and 1b (charged). If the Injector is charged and the Handles are closed, proceed to **Section 2) Fill Syringe**. DO NOT open the Handles when the Injector is charged.

figure 1a

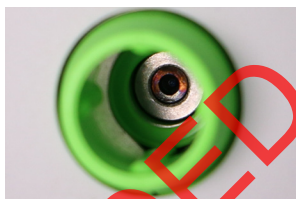
figure 1b

- b) When the Injector is in the uncharged state, as shown in figure 1a, and the Handles are closed, open the Handles to prepare the Injector. Press the rear Injector Latch Button to allow the Injector Handles to open freely, as shown in figure 2.

figure 2

NOTE: If the Injector is in the charged state, as shown in figure 1b, and the Handles are partially open, as shown in figure 3, close the Handles and proceed to **Section 2) Fill Syringe**. DO NOT touch internal components or attempt to fully open the Handles when the Injector is charged.

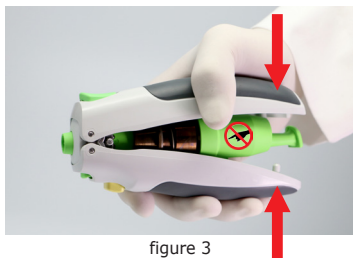
figure 3

- c) With the Injector Handles fully open, as shown in figure 4, squeeze the Handles until they close and are locked in place, as shown in figure 5. Once the Injector Handles are closed, ensure the Injector Hammer is in the charged state, as shown in figure 1b.

figure 4

figure 5

NOTE: If the Injector continuously fails to reset, refer to the contact information found at the end of this document and contact a PharmaJet representative.

2) Fill Syringe

- a) Remove protective cap from the primary vial packaging and prepare the top following aseptic techniques.
- b) Aseptically remove the Filling Adapter from its packaging.
- c) Aseptically remove the Syringe from its packaging. Remove Protective Cap.
- d) Insert Syringe into Filling Adapter by employing a clockwise turn, as shown in figure 6. Slide Protective sleeve off of the Filling Adapter spike.
- e) Insert the Filling Adapter spike into the vial with a twisting motion until it is fully inserted, as shown in figure 7.
- Note:** Adding air to a vial can minimize vacuum and assist filling.
- f) Invert the assembly.
- g) Pull down Plunger to draw fluid into Syringe. Remove bubbles as required. Complete filling by aligning red rubber ring with the fill line, as shown in figure 8.
- h) Snap off end of Plunger protruding from Syringe, as shown in figure 9.

figure 6

figure 7

figure 8

figure 9

Note: The Filling Adapter is for single vial use. It may be used to continuously fill several Syringes from the same multi-dose vial for up to 48 hours. Up to five vial punctures with a new Filling Adapter can be performed on a multi-dose vial. Cap the Filling Adapter with a new unused Syringe between fillings.

Note: Prior to vial storage, remove and dispose of the Filling Adapter in accordance with local requirements.

3) Load Injector

- a) Insert the Syringe and Vial Assembly into the Inner Core Housing of the charged Injector until it “clicks” into place, as shown in figure 10.

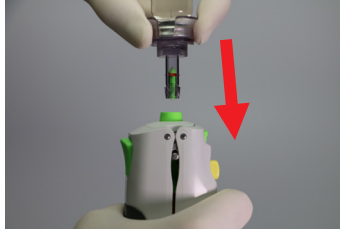
figure 10

- b) Invert and remove the Filling Adapter & vial from Syringe by employing a counterclockwise turn.

4) Give Injection

NOTE: If needed, remove protective Syringe Cap.

- a) With forefinger, depress and hold the Injector Safety Button. This allows the Inner Core Housing to retract freely, as shown in figure 11.

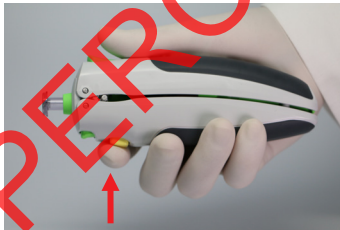
figure 11

- b) While holding the Safety Button, press the Syringe perpendicular to the injection site, as shown in figure 12. Once the Inner Core Housing is retracted, it is not necessary to hold the Safety Button as long as you keep pressure on the patient and the Inner Core Housing remains retracted.

figure 12

c) To inject, depress Activation Button with thumb. A “click” will be heard upon injection.

d) Verify Syringe contents were delivered by observing that the Plunger has advanced completely to the Syringe tip.

NOTE: If the Inner Core Housing returns to the original extended position and becomes locked, simply press the Safety Button again (See Figure 11) and start injection operation from Section 4.

5) Eject the Syringe

a) Remove Syringe from the Injector by pressing the Eject button, as shown in figure 13.

figure 13

b) Dispose of Syringe in accordance with local requirements.

c) For storage, keep the Injector in the closed configuration.

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